


## LIASYS 330 DURABLE CUVETTE SYSTEM

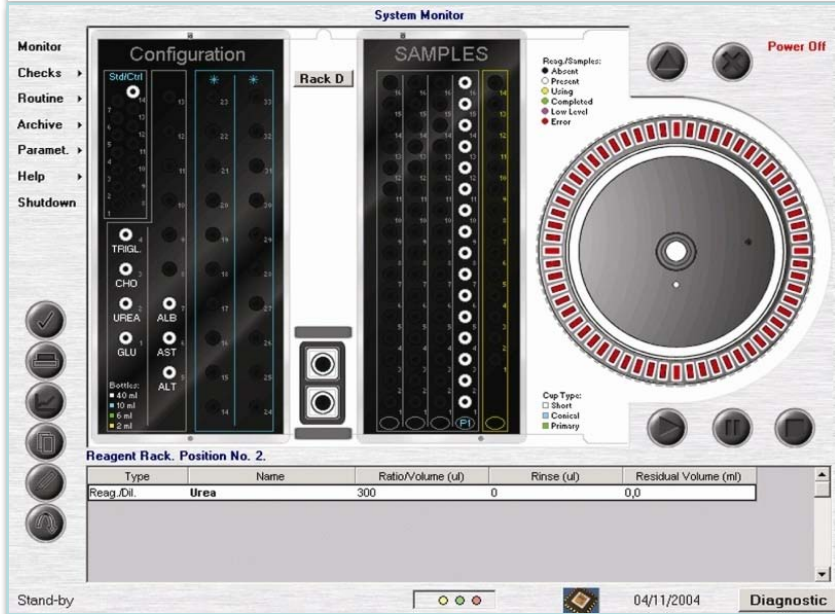

60 optical grade cuvettes individually set on in a high capacity reaction wheel, are continuously tested for cleanliness to eliminate carryovers.

WBL tests (Water Blank Level) run by the on-board Quality Assurance Program monitor each cuvette and warn if any of the segments includes a cuvette that does not meet the minimum quality requirements; if needed, each cuvette segment can be easily replaced by the system operator.

The intuitive and simple graphical interface provides immediate system and operation feed-back; like durable cuvette status, test results time toward tests or work-list completion.



L 330 Reaction wheel with durable cuvettes

DURABLE OPTICAL CUVETTES - QC MONITORED AND USER CHANGEABLE